



COLINK[®]

Standard & Adjustable

Fixation button
with loop



Teknimed

COLINK[®] Standard

Fixation button with loop

(Class IIb)

COLINK[®] Standard fixation device offers one of the strongest soft tissue femoral fixation currently available.

MECHANICAL AVANTAGES:

- **HIGH TENSILE STRENGTH (1500N)¹**
- **CONTINUOUS BRAIDED LOOP**

HANDLING AVANTAGES:

- **TRULY ENDOSCOPIC PROCEDURE**
- **ACCOMODATES VARIOUS GRAFT LENGTHS**
- **ELIMINATES THE NEED FOR KNOT TYING**
- **ENSURES CORRECT SEATING OF DEVICE**

Ideal for primary or auxilliary fixation ACL or PCL reconstruction techniques. Preloaded with UHMWPE suture (white) and flipping suture (striped) for added procedure efficiency.

Can be used in Single-Bundle soft tissue fixation; Double-Bundle soft tissue fixation; Bone-Tendon-Bone fixation.

Indications


Use in the fixation of bone and soft tissue in orthopedic procedures requiring ligament or tendon reconstruction.

Composition

Loop: UHMWPE

Sutures: HS Fiber[®]

Plate: Titanium (Ti-6Al-4V)

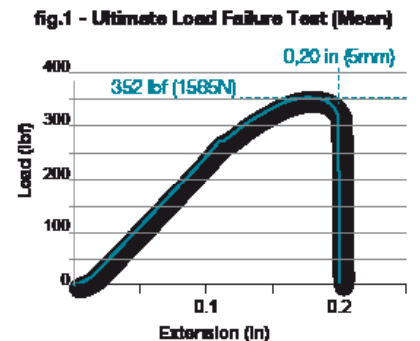

Endoscopic Reamer Ø 4.7

AVAILABLES SIZES

Standard 12 mm	Standard 35 mm
Standard 15 mm	Standard 40 mm
Standard 20 mm	Standard 45 mm
Standard 25 mm	Standard 50 mm
Standard 30 mm	Standard 55 mm

1, Data on file at Teknimed.

COLINK[®] Adjustable

Fixation button with loop

(Class IIb)

COLINK[®] Adjustable fixation device offers leading strength and incorporating a unique cradle design to help protect the graft during loop reduction.

MECHANICAL AVANTAGES:

- **HIGH TENSILE STRENGTH (900N)**
- **UNIQUE 3-POINT LOCKING DESIGN & BRAID FOR EASY PLACEMENT AND TO REDUCE SUTURE CREEP**
- **GRAFT SUPPORT FRAME WITH LARGE CONTACT AREA**
- **ENSURES PROPER REDUCTION & COUNTER TENSIONING**

HANDLING AVANTAGES:

- **ADJUSTABLE LOOP - ONE SIZE FITS ALL (15 to 60mm)**
- **TWO-HANDED ADJUSTMENT**
- **CROSSOVER PATTERN TO INCREASE VISIBILITY**

Indications

Use in the fixation of bone and soft tissue in orthopedic procedures requiring ligament or tendon reconstruction.

Composition

Loop: UHMWPE
Sutures: HS Fiber[®]
Plate: Titanium (Ti-6Al-4V)



Graft Prep Station



COLINK®

Storage

COLINK® is supplied sterile, ready to use in operating rooms. It must be kept in its original unopened packaging, in a dry, clean place away from light and at room temperature.



Sku

COLINK® Standard	
SIZE	SKU
Colink® standard 12 mm	TEK-OBCL-12
Colink® standard 15 mm	TEK-OBCL-15
Colink® standard 20 mm	TEK-OBCL-20
Colink® standard 25 mm	TEK-OBCL-25
Colink® standard 30 mm	TEK-OBCL-30
Colink® standard 35 mm	TEK-OBCL-35
Colink® standard 40 mm	TEK-OBCL-40
Colink® standard 45 mm	TEK-OBCL-45
Colink® standard 50 mm	TEK-OBCL-50
Colink® standard 55 mm	TEK-OBCL-55
Colink® standard 60 mm	TEK-OBCL-60

COLINK® Adjustable
TEK-OBAL-80

COLINK® Graft prep station
BAK-7069

COLINK® Endoscopic Reamer	
Ø 4.7mm	BAK-7118


MD Class IIb

CE 0459



LEGAL MANUFACTURER:
Riverpoint Medical
 825 NE 25th Avenue
 97232 Portland
 United States

Distributed by

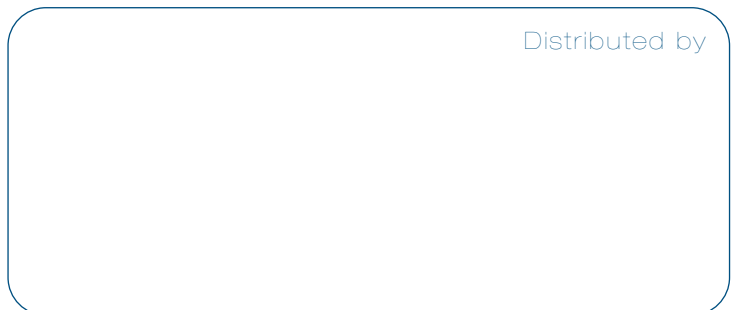

Sales and engineering departments:

Z.I. de Montredon - 12, rue Apollo
 31240 L'UNION (France)
 Tel (33) 5 34 25 10 60
 Fax (33) 5 34 25 27 39

Administration and supply-chain offices:

Z.I de la Herray
 65500 VIC en BIGORRE (France)
 Tel (33) 5 62 96 88 38
 Fax (33) 5 62 96 28 72

